





Innovation in Gynecology

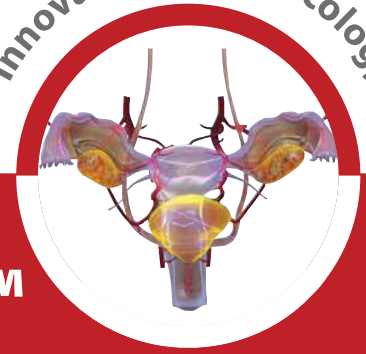


LiNA McCartney Tube™


*Improving the outcome of
Minimally Invasive Gynaecology*

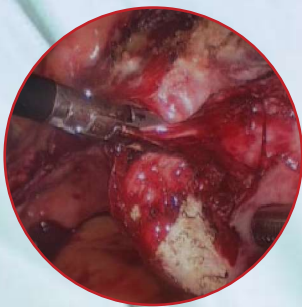


Identifies where to safely divide the uterine vessels

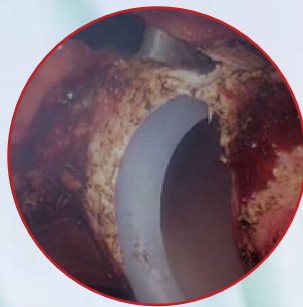

Provides a clear template for safe and easy colpotomy

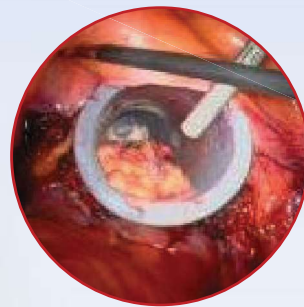

Maintains pneumoperitoneum and provides an easy exit for specimens

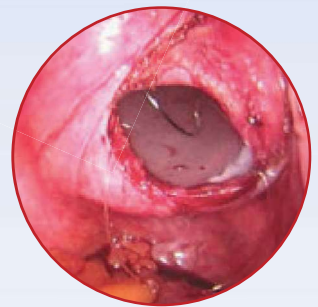

Delivers a convenient entry and exit for needle and suture and presents the vaginal edges for laparoscopic closure of the vaginal vault



LiNA McCartney Tube™

The LiNA McCartney Tube™ offers you a unique single use device, sterile and ready to use.



safe

- Maintains pneumoperitoneum after colpotomy
- Protection of the ureters from trauma during surgery
- Identifies where to safely seal and divide the uterine vessels
- Provides a clear template for bladder dissection and colpotomy
- Improved visualisation designed for a safer procedure

functional

- **Lid with 5mm & 10mm valves**
Prevents loss of pneumoperitoneum and enables insertion of instruments.
- **Fast and simple laparoscopic suturing technique of the vaginal vault**
Delivers a convenient entry and exit for needle and suture and presents the vaginal edges for laparoscopic closure.
- **Handle to control tube positioning during procedures**



versatile

- **Easy removal of specimen, e.g. the uterus, lymph nodes and ovarian cysts**
- **Tip withstands both electrosurgical and ultrasonic energy**
Soft curvature designed for atraumatic insertion and anatomic fit.




Ordering Information:


LiNA McCartney Tube™  2797

35-TV	LiNA McCartney Tube™ - 35mm diameter, sterile, single use	Box of 5
45-TV	LiNA McCartney Tube™ - 45mm diameter, sterile, single use	Box of 5

© 2019 LiNA Medical ApS. All rights reserved.

EN ISO 13485:2016
Quality Management System

 MANUFACTURER: The OR Company Pty Ltd 1/32 Silwood Rise Carrum Downs VIC 3201 Australia

 OBELIS S.A., BD GENERAL WAHIS 53, B-1030 BRUSSELS, BELGIUM, MAIL@OBELIS.NET

■ Exclusive Distributor for Designated Countries: LiNA Medical International Operations AG, Platz 3, 6038 Root D4, Switzerland.